


눈감기->클릭->고급클릭

FrameRate: 61.39640

루지통

Basic Mouse Control

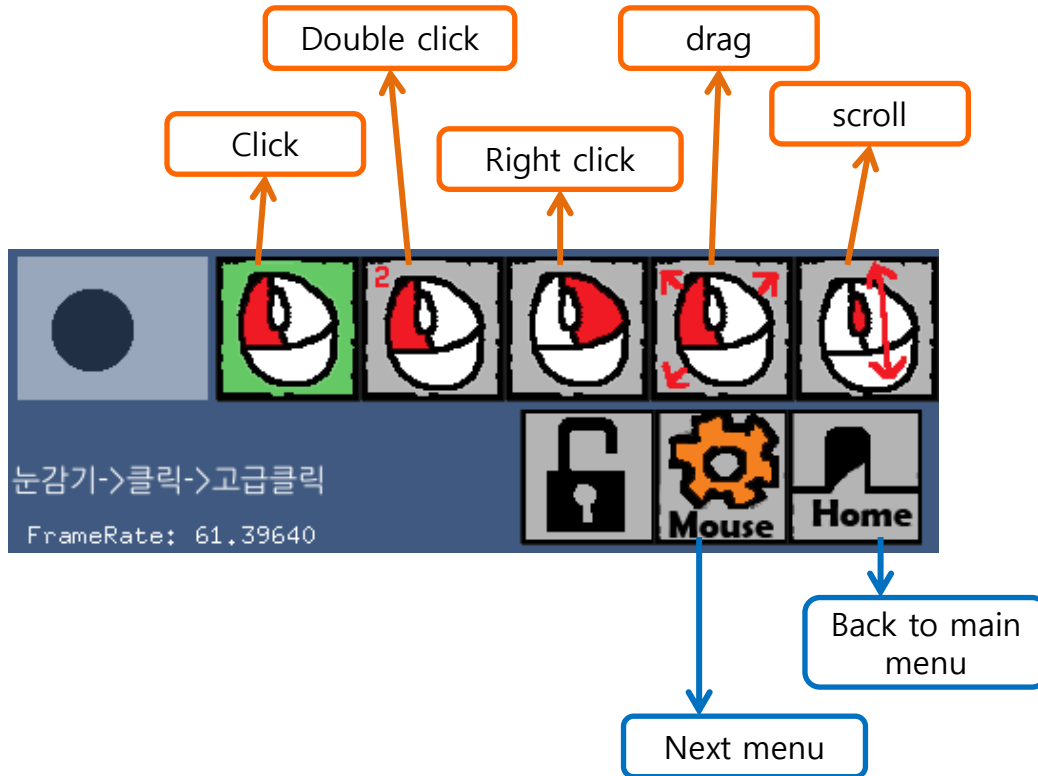


Click the click action you want to perform
If you click twice, the setting will be fixed.





<Action state fixed to right click>

Changing Mode

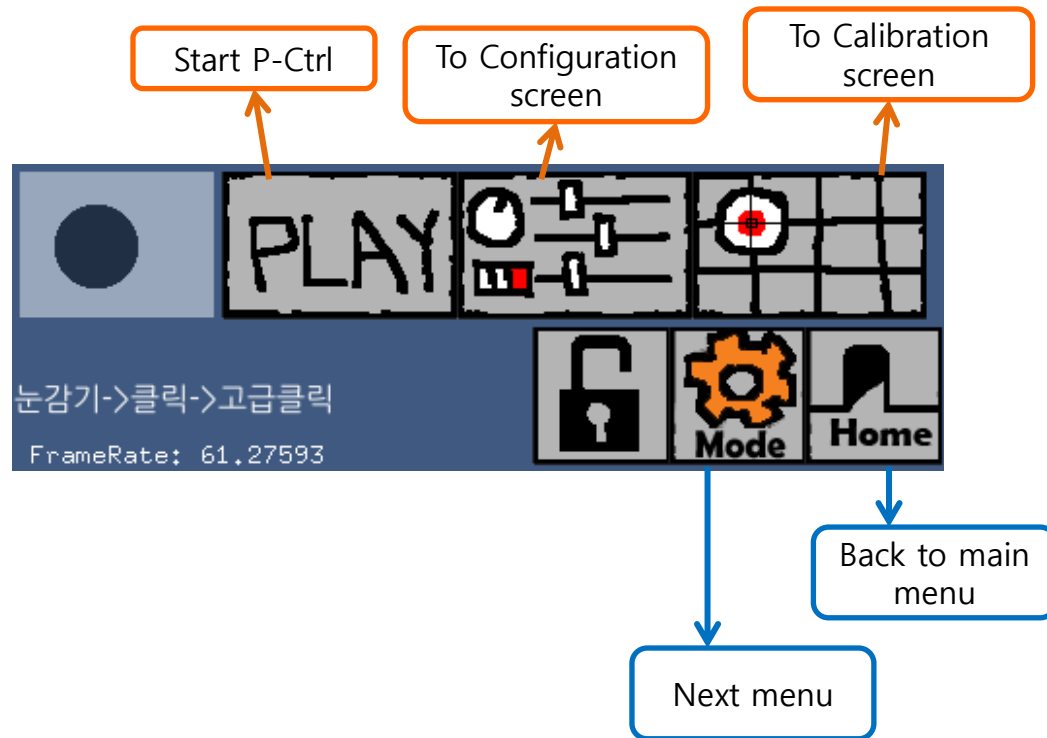

Hotkeys

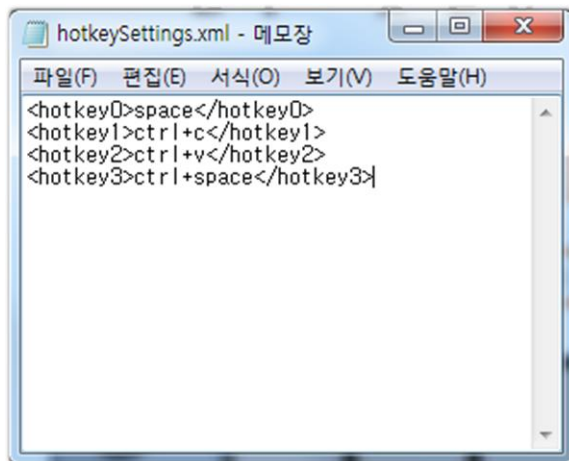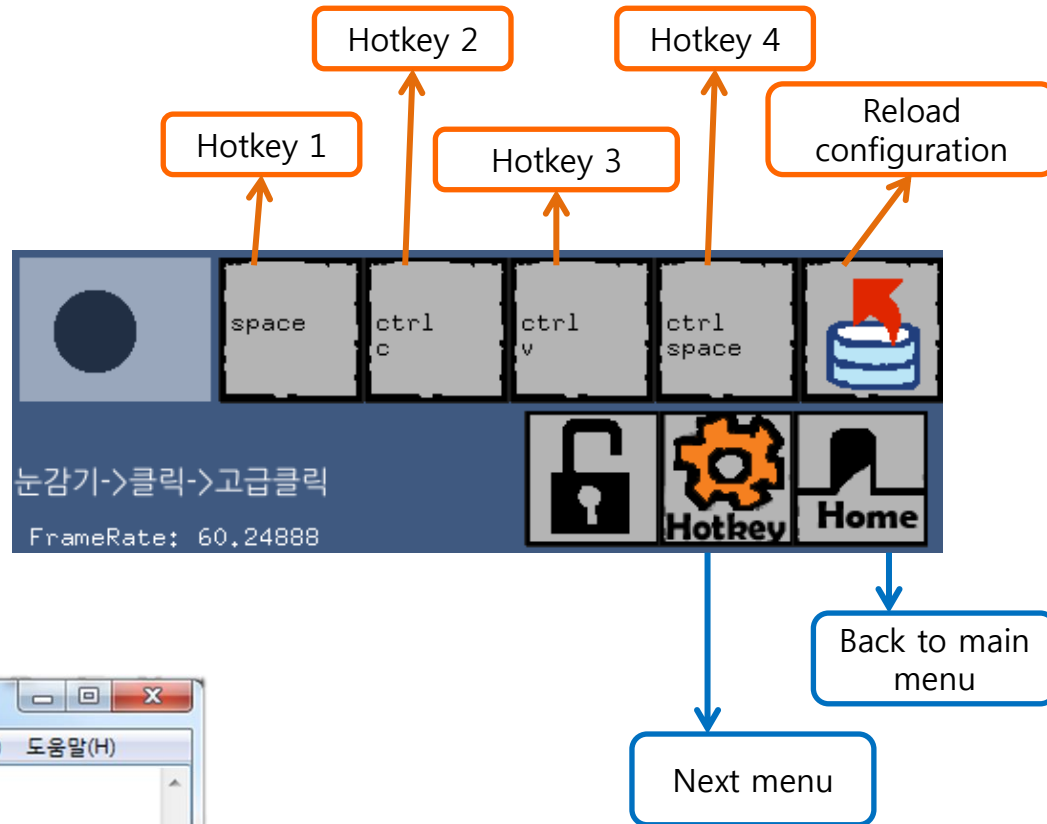

In (eyeCan program folder)\data\settings folder, edit the hotkeySettings.xml file to change the hotkey settings.

Lock Mode



Lock




Lock to block any click events to the PC.
This may be useful when watching videos

Advanced Click

How to
: perform long click (long blink)
when using mouse

Description

: A box will appear on the screen.
You can move the mouse cursor in
the box by looking at 4-directions
around the box. You may want to
use this to perform click actions
that requires greater accuracy.
Long-click to cancel or go to Error
compensation mode.

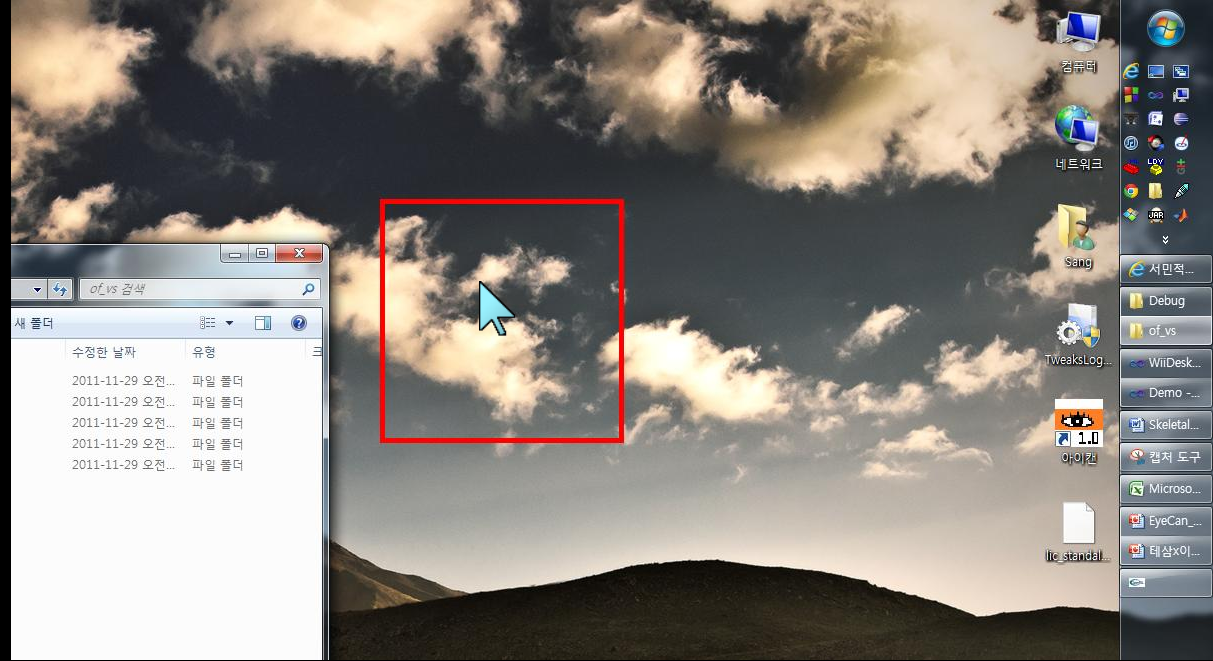

Advanced Uses

Error Compensation

How to
: perform long click (long blink) when
using Advanced Click mode

Description

: When your calibration has gone
wrong, you can make mouse cursor
position offset to recover this
calibration-off problem. Look at the
center of the box for a while (for a
few seconds) then the program will
automatically compute the
appropriate offset.
Long-click again to cancel.

